




THUASNE

Thuasne USA 3D Scanning

Scan for your custom braces
in under 60 seconds



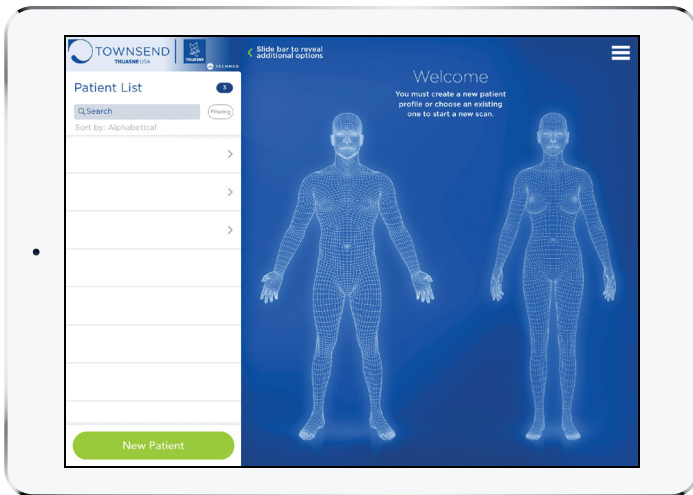

Thuasne USA 3D Scanning

- TechMed 3d technology (Industry Standard)
- Precise 3D scans in seconds
- Quicker turn-around times
- Consistent Quality
- Thousands of Optimal Scans received
- Lightweight & easy to use
- No charge for file conversion (Thuasne USA Only)
- All extremities
- Files Stored in perpetuity (Forever!!!)

Compatibility

- Apple products only at this time

**recommended use of I pad for easier process & success*

Thuasne USA Models Available from Scan

Ligament Knee	OA Knee	Ankle/ Foot & Elbow	Specialty KO	KAFO
Rebel Line (Lite, Std, Pro)	Rebel Reliever	Premier Ankle	Full Shell Trigger Lock	Premier KAFO
Premier	Premier Reliever	Premier AFO	Premier Trigger Lock	Full Shell KAFO
Air Townsend	Premier Reliever 1	<i>SpryStep® Custom Coming 2019</i>	Full Shell Polio	Hybrid KAFO
Air Townsend	Lite Air Townsend	<i>SpryStep® Vector Coming 2019</i>	Premier Polio	<i>Dynamic KAFO Coming 2019</i>
	Air Townsend Reliever lite	Full Shell Elbow	Full Shell BK	
	UniReliever			
	Reliever One			

Optimal Scanning Positions

Hinged AFO – 90° @ Ankle, 90° @ Knee, 90° @ hip

Non-Hinged AFO – 90° TKA/ Weight Line & 7° Tibial Inclination

KO – Leg fully extended, foot dorsi flexed

Tips & Info

- Purchase App Thuasne 3D in the App store on your apple device (it's FREE)
- All files are sent to escan@thuasneusa.com
- You must designate this your "email recipient" in the settings within the App
- "Scan Box Placement Strategy" within App should be set to "Automatic" for KO & AFO purposes
- You must elevate the leg or angle of the footplate while taking scan (Thuasne USA Scanning Platform)
- "Scan Box Adjustment" option allows you to capture a scan from a vertical or horizontal view
- Attach Detailed Notes & Pictures
- If you mess up, don't be afraid to start over and try again
- Send us multiple scans if you like, we will pick the best one (no extra charge)
- Scanning Maximum Circumference Capability is 35"

Pricing

- Structure Scanner \$349 USD
- Structure Scanner w/ Bracket \$379 USD
- Bracket (Model specific) \$50 USD



To purchase or shop scanner & accessories, please go to:

<http://store.structure.io/buy/structure-sensor>

Thuasne USA Scanning Platform

- \$275 CAD



5545 St. Jacques O. Montreal, QC H4A 2E3
 Local: 514. 369. 3311
 Toll Free / Sans Frais : 1.877.222. 3311
 Fax Toll Free / Telec Sans Frais : 1.877.527.1911
www.bmiortho.com