

Premium soft braces The optimum in functional excellence, comfort and design.



The exclusive Canadian distributor for Thuasne USA



Clear structures, innovative solutions The new product range from THUASNE THÄMERT.

Successful partnerships between trusted medical stockists and manufacturers will increasingly be dictated by the ability to support the active selling process with a range of products that perfectly matches requirements and is segmented accordingly.

In addition to supplying an affordable core range, offering enhanced products for customers demanding the highest standards in terms of design and comfort is becoming increasingly important.

With the new **PR0**MASTER range of premium soft braces we are consciously taking one step further towards offering our customers a clearly structured product concept that is perfectly tailored to market requirements.

PROMASTER sets new standards in terms of functional excellence, design and quality. As the flagship of our new design language, it clearly mirrors the changes taking place in our entire company.

Together with leading industrial designers, our in-house development team has worked closely with our international Group to make sure this new range of products provides both active athletes and recuperating patients with modern, comfortable and effective care.

We hope you enjoy reading this brochure. Your THUASNE THÄMERT

Keep moving!

PROMASTER premium soft braces Outstanding design. Innovative details.

What makes the new flagship product developed by THUASNE THÄMERT so special? Its professional design, superb quality and innovative technology that sets new standards.

Design

Together with THUASNE THÄMERT experts, well-known industrial designers have developed a subtle, contemporary design that is superbly functional and will appeal to both trend-oriented patients and active athletes alike.

FlexAir pads

Pads with a new structure have been developed specially for the **PR0**MASTER range: they adjust snugly to the shape of the body and are point elastic. Trigonal elements spaced apart from each other enhance the flexible bending properties of the FlexAir pads, allowing them to adjust perfectly to each body movement. A ventilation system transports moisture away from the skin, thus making the soft braces even more comfortable.

Proprioception

While shaping and designing the surface texture of the pads, particular attention was paid to stimulating the muscle-tendon junctions. For example, nubs in sensitive joint areas reinforce proprioceptive perception.

- Breathable knit fabric with 3-D texture and latex-free elastane guarantee a pleasantly firm fit.
- Sturdy grips make them easy to put on.
- Gradual change in compression towards the edges
 The knit pattern has been optimised to ensure that compression decreases from the centre of the joint to the edge of the brace.
 The resulting soft, wave-shaped change in compression guarantees a perfect fit and stops the braces from pinching.

Genu PROMASTER Stimulates the muscles. Relieves pressure points.

The soft knee brace features adaptive pads that specifically stimulate proprioceptive areas to correct posture. The knee is kept stable, the patella supported even more reliably, and pain is alleviated.

- Y pad for precision patella hold
- Integrated nubs stimulate the thigh muscles and accelerate the healing process
- FlexAir structure on the pad ensures that the brace fits perfectly around the patella, and that air can circulate





Y pad

- Patella support
- Surrounding tissue is massaged, extending right up to the thighs
- The distal area of the pad is cut out to relieve pressure on the edge of the shinbone



Comfortable, firm fit

- Gradual change in compression prevents pinching
- Breathable knit fabric with 3-D texture for comfort




Grips on both sides

- Make it easier to apply the soft brace
- Patent pending



Welded spiral springs

- Enhanced stabilisation of the knee joint
- Gently support knee movement, without obstructing it
- Stop the brace from slipping

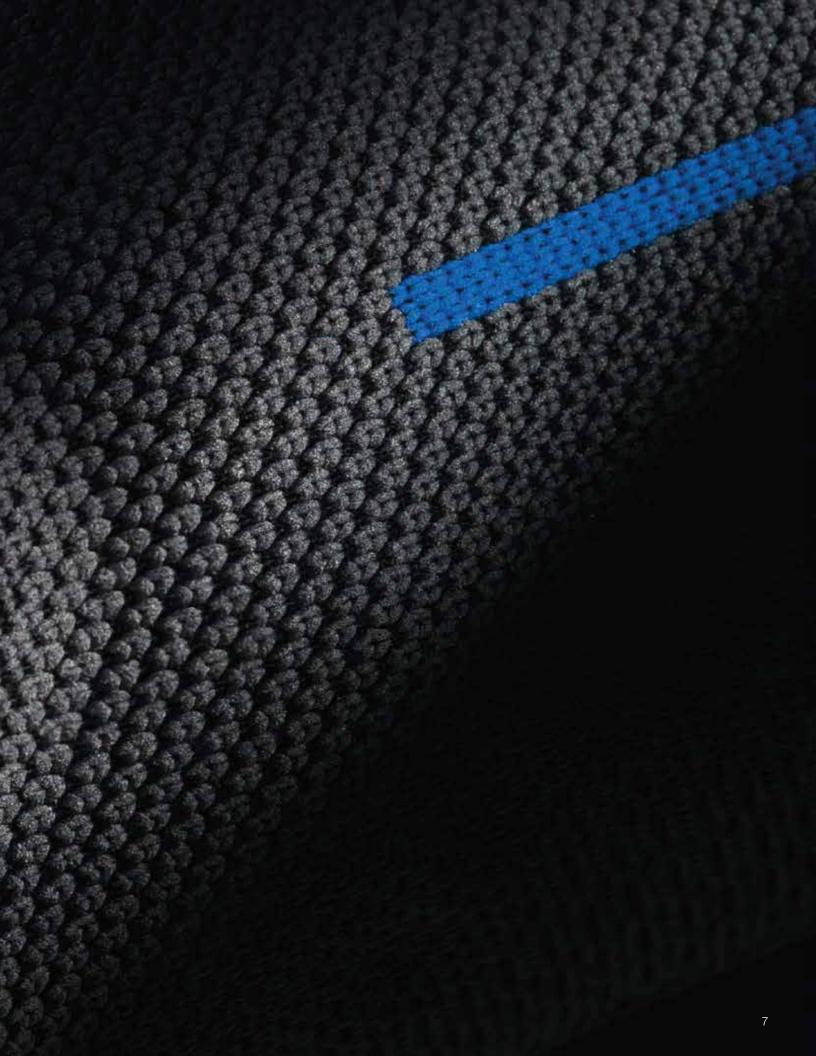

Manu PRO MASTER Superb mobility. Outstandingly comfortable.

Manu PR0 MASTER relieves and supports the wrist following injuries, or when patients are suffering from carpal tunnel syndrome or wrist arthrosis. The point-elastic FlexAir structure on the pad channels pressure away from nerves and blood vessels. The soft-knit thumb joint ensures maximum comfort and mobility.

- Pain following post-operative or post-traumatic irritation is relieved quickly and lastingly
- Supports and stabilises the joint
- · Adapts to any hand movement without slipping
- The compression effect of the knit fabric with 3-D texture improves blood circulation in the injured wrist



Soft-knit joint

• Maximum comfort at the junction between thumb and index finger



Elastic compression strap and plantar splint

• Enhanced stabilisation of the wrist



Grip on the inside arm

• Helps when putting the soft brace on

FlexAir pad

- Relieves pressure on the nerves
- Alleviates pain
- Adjusts snugly to the shape of the body

Malleo PRO MASTER Firm hold. All-round support.

The anatomically shaped **Malleo PRO** MASTER supports the ankle on both sides. Even compression, supported by point-elastic pads, reduces oedema and pain. Soft, highly flexible knit fabric for enhanced comfort.

- Air-permeable, point-elastic silicone pads reduce oedema
- Breathable fabric
- Decreased compression at the edges of the soft brace prevent pinching
- Available for left or right ankles







Anatomically shaped pads

- Effectively support the ankle
- FlexAir structure guarantees the circulation of air



Patented pronation strap on Malleo PRO MASTER Aktiv

- Enhanced forefoot stability
- Additional protection against the foot buckling inwards
- Thus preventing resulting additional stress on the ankle ligaments and tendons



Epi PROMASTER Wonderfully comfortable. Stimulatingly effective.

Epi PROMASTER reduces and eliminates tension and pain in the elbow muscle insertions. When the elbow is moved, the point-elastic FlexAir structure of the pad activates the musculature and thus optimises the healing process.

- Alleviates pain
- Flexible, point-elastic pads stimulate the musculature
- Soft, 3-D knit fabric adapts to each movement of the body
- Promotes healing and mobility





Anatomical fit

- Supports the joint areaPrevents slipping



Soft knit fabric

• Maximum comfort in the elbow crease



FlexAir pads

• Relieve pressure on the tendon insertions in the forearm muscles





Genu OA

Supportive and stabilising orthosis for gonarthrosis patients. Slim yet powerful.



Indications	Gonarthro	Gonarthrosis; to partially support the medial or lateral compartment						
How it works		Uses the 3-point principle to support one compartment and physiologically stabilises the knee joint with adjustable restriction of flexion/extension.						
Colour	Anthracite	5						
Size	Size XS S M L XL 2XL 3XL 4XL medial rig	SizeThigh circumference inchCalf circumference cmXS $38 - 41$ $15 - 16$ $28 - 31$ $11 - 12 1/4$ S $41 - 44$ $16 - 17 1/4$ $31 - 34$ $12 1/4 - 13 1/2$ M $44 - 47$ $17 1/4 - 18 1/2$ $34 - 37$ $13 1/4 - 14 2/3$ L $47 - 50$ $18 1/2 - 19 3/4$ $37 - 40$ $14 2/3 - 15 3/4$ XL $50 - 53$ $19 3/4 - 21$ $40 - 43$ $15 3/4 - 17$ 2XL $53 - 57$ $21 - 22$ $43 - 46$ $17 - 18$ 3XL $57 - 61$ $22 1/2 - 24$ $46 - 49$ $18 - 19 1/3$					•	
ProdCode	35 047							

Slim yet powerful



Innovative, effective, comfortable:


Silicone threads in the top knitted edge

lessen the threat of the orthosis slipping.



Both medial and lateral arthrosis can be treated

e.g. medial right/later left or medial left/lateral right.

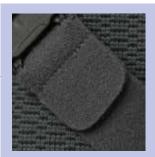

Individually adjustable straps on the front and back of the orthosis

for personalised adjustment to ensure the orthosis fits firmly.



Cross-over straps

for dedicated support and optimised pressure distribution.



Flat joint developed specially for the orthosis

fits close to the leg without being conspicuous.

Flexion and extension settings

Flexion settings: 20°, 40° Extension settings: 0°, 10°



Individually adjustable angles

for an optimised fit.



Magnetic fasteners at the front

make the straps easy to open and close.

Genu PROMASTER Soft knee brace with reinforced proprioceptive effect

Particular features	Indications	How it works
 Y-shaped patella pad 3-D knit fabric covering the patella and in the hollow (highly flexible) Silicone nubs at the sides to prevent slipping Latex-free Springs integrated at the sides 	 Post-traumatic, post-operative and chronic irritation Tendomyopathy and ligament insertion degeneration Arthrosis and arthritis Minor joint instability Patella pathologies 	Provides knee compression and support

Manu PROMASTER Soft wrist brace with reinforced proprioceptive effect

	Particular features	Indications	How it works
2	 Elasticated knit brace 3-D knit fabric with pad cutout over the condyle Stabilising splint Compression strap Latex-free 	 Tendovaginitis Post-traumatic, post-operative and chronic irritation Sprains Joint effusions and swellings caused by arthritis and arthrosis 	Provides wrist compression and stabilisation

Malleo PRO MASTER soft ankle brace with reinforced proprioceptive effect

Particular features	Indications	How it works
 Elasticated knit brace 3-D knit fabric over the ankle bones Malleolus pads Latex-free 	 Post-traumatic, post-operative and chronic irritation Sprains Joint effusions and swellings caused by arthritis and arthrosis Minor joint instability 	Provides ankle compression and support

Malleo PRO MASTER Aktiv soft ankle brace with reinforced proprioceptive effect and pronation strap

Particular features	Indications	How it works
 Elastic knit brace with integrated pronation strap 3-D knit fabric over the ankle bones Malleolus pads Latex-free Available from autumn 2014 onwards 	 Post-traumatic, post-operative and chronic irritation Sprains Joint effusions and swellings caused by arthritis and arthrosis Minor joint instability 	Provides ankle compression and support

Epi PROMASTER soft elbow brace with reinforced proprioceptive effect

Particular features	Indications	How it works
Elasticated knit brace	• Epicondylitis	Provides elbow
 Pad cutout over the condyle 	 Post-traumatic, post-operative and 	compression and
Epicondyle pads	chronic irritation	stabilisation
Latex-free	• Sprains	
	 Joint effusions and swellings caused by 	
	arthritis and arthrosis	

Reversa PRO MASTER soft lumbar brace with reinforced proprioceptive effect

Particular features	Indications	How it works
 Elasticated knit brace 2 dorsal and 2 lateral reinforcement rods Pocket for lumbar pad Additional grip makes it easier to put on Latex-free Available with silicone or Flex pad 	 Post-traumatic, post-operative and chronic irritation Lower back pain with or without acute or persistent sciatica Lumbar spine support following surgery Connective tissue diseases Muscle fibre pain 	Supports and stabilises the spine

Colour	Version	Dimens	sions		Art. no.	HMV. no.	MiGel no.
Dark grey/ Blue	Symmetric	Size XXS S M H L L+ XL+ XL+ XL+ XXL+	$\begin{array}{c} \text{Thigh circumference} \\ 35-38 \ \text{cm} \\ 38-41 \ \text{cm} \\ 41-44 \ \text{cm} \\ 44-47 \ \text{cm} \\ 47-50 \ \text{cm} \\ 47-50 \ \text{cm} \\ 50-53 \ \text{cm} \\ 50-53 \ \text{cm} \\ 53-56 \ \text{cm} \\ 53-56 \ \text{cm} \\ 56-59 \ \text{cm} \end{array}$	Calf circumference 25-28 cm 28-31 cm 31-34 cm 34-37 cm 37-40 cm 37-40 cm 40-43 cm 40-43 cm 43-46 cm	35 040	05.04.01.0118	05.04.02.00.1

Colour	Version	Dimensions	Art. no.	HMV. no.	MiGel no.
Dark grey/ Blue	Right or left	Size Wrist circumference XS 14–15 cm S 15–16 cm M 16–17 cm L 17–18 cm XL 18–19 cm XXL 19–20 cm	35 070	05.07.02.0143	05.07.03.00.1

Colour	Version	Dimensions	Art. no.	HMV. no.	MiGel no.
Dark grey/ Blue	Right or left	Size Ankle circumference XS 17–19 cm S 19–21 cm M 21–23 cm L 23–25 cm XL 25–27 cm XXL 27–29 cm	35 020	05.02.01.0075	05.02.02.00.1

Colour	Version	Dimensions	Art. no.	HMV. no.	MiGel no.
Dark grey/ Blue	Right or left	Size Ankle circumference XS 17–19 cm S 19–21 cm M 21–23 cm L 23–25 cm XL 25–27 cm XXL 27–29 cm	35 022	05.02.01.2XXX (application filed)	05.02.02.00.1

Colour	Version	Dimensions	Art. no.	HMV. no.	MiGel no.
Dark grey/ Blue	Symmetric	Size Forearm circumference XS 21-23 cm S 23-25 cm M 25-27 cm L 27-29 cm XL 29-31 cm XXL 31-34 cm	35 080	05.08.01.1053	05.08.02.00.1

Colour	Version	Dimensions	Art. no.	HMV. no.	MiGel no.
Dark grey/ Blue	Symmetric	Size Body circumference XS 70- 80 cm S 80- 90 cm M 90-100 cm 100-110 cm L 100-120 cm XL XXL 120-130 cm 120-130 cm	35 140 Straight 35 145 Fitted	23.14.03.2078	05.14.04.00.1





5545 St. Jacques O. Montreal, QC H4A 2E3 O: 514.369.3311 Toll Free / Sans Frais : 1.877.222.3311 Fax Toll Free / Telec Sans Frais : 1.877.527.3311 www.bmiortho.com